



## HEALBOND MP

Universal self-etching bonding agent

**TEGDMA, Bis GMA, HEMA-free**

### Features


- Multi-purpose (direct and indirect restorations) and universal : excellent shear bond strength on tooth structures, ceramic, metal and Zirconia.
- Good wetting properties to be used on slightly moist dentin surfaces (wet-bonding-technique).
- Contains Benzalkonium Chloride to eliminate residual bacteria (no post-operative sensitivities).
- Releases fluoride.

### Indications

- Light-curing composites on tooth structures.
- Light-curing composites on metal (replacement of broken ceramic chips in mouth).
- Dual-curing composite (core build-up), when mixed with Healbond A activator.

### In conjunction with Innocem :

- Ceramic veneers, root canal posts and Maryland bridges.

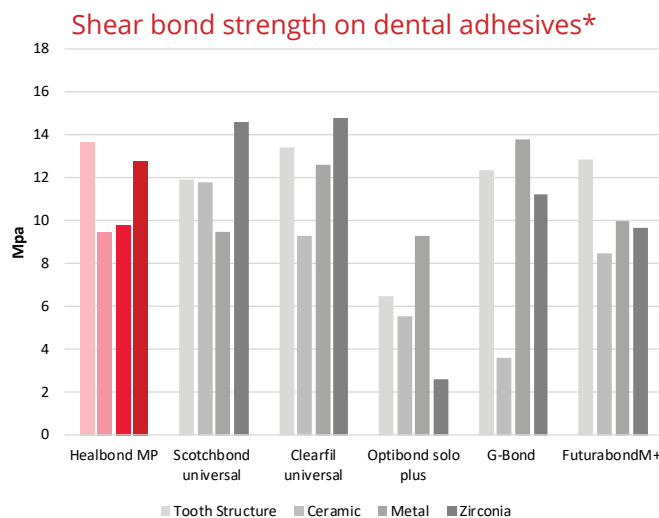

### References & presentations

BMP-5 5ml bottle.

BMPA-5 Activator for dual-curing use. 5ml bottle.

MD : Class IIa certified by Tüv Rheinland.

\* Tests conducted by G-Pharma. Shear bond strength : tests carried out following the ISO 11405 norm. Conversion rate : the degree of conversion was evaluated using a FTIR spectrometer with an attenuated total reflectance (ATR) accessory. The polymerization has been conducted by an Elsodent lamp: 9mm of diameter and 600mW/cm<sup>2</sup> power at room.



**Elsodent Laboratory**

35 Avenue du Gros Chêne - BP 50072 - Herblay 95612 Cergy Pontoise Cedex - France

Phone: +33 (0)1 30 37 75 75

elsodent@elsodent.com - www.elsodent.com